MACHEREY-NAGEL NANOCOLOR® VIS II and NANOCOLOR® UV/vis II







CIE-conform color measurements and analyses


Integrated inspection equipment monitoring tools



Reach your goals smartly

Clearly arranged icons instead of complicated lists

With the brand-new spectrophotometers *NANOCOLOR® VIS II* and *NANOCOLOR® UV/VIS II*, MACHEREY-NAGEL revolutionizes your daily laborato- ry work, combining premium high-tech instruments with outstanding usability. Work with these innovative photometers as with your smartphone or tablet and enjoy the intuitive, icon-based menu. The clearly arranged and high-resolution touch screen display makes your daily measurement routine a real pleasure.

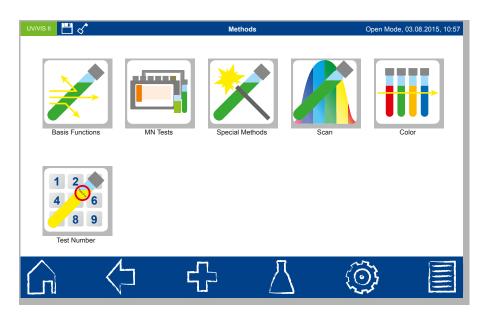

The __icon allows a quick access to all programmed methods.



User guidance

- Smart, icon-based menu guidance
- Easy to use without time-consuming training
- Intuitive operation for fast results
- → Revolutionary user experience



Touch screen

- Projective-capacitive touch screen
- 10inchHDtouchscreendisplayforcomfortablelabwork
- Reliable operation even with gloves
- Innovative usage concept



Makes your colleagues envious

Powerful technology behind the scenes

The brand-new NANOCOLOR® spectrophotometers impress with high-class technology and optics. The spectral bandwidth of the NANOCOLOR® UV/v/s II has been reduced to less than 2 nm and allows high-precision measurements. The wavelength resolution of 0.1 nm enables a new level of precision, for example for the monitoring of high-resolution scans. Internal stray light is reduced to < 0.05 %, resulting in extremely accurate measuring signals. The optical setup and the clever technology of both instruments enable measurements without protective cover, a big advantage for seamless lab processes.





Universal cuvette holder without cuvette slot cover

- Cuvette slot for 10, 16, 20, 40, 50 mm cuvettes
- Measurements without cuvette adapter and slot cover
- Insensitive to external light
- Automatic detection of cuvette size
- Seamless laboratory processes



Sample detection and sample information

- Time saving sample allocation with a handheld scanner
- Quick and easy administration of all important sample information
- Future-proof data management



Optical setup

- Precise optics with reference detector technology
- High-resolutionscans
- Minimized aberrations due to holographic gratings
- Strong technology for more precise measurements



2D barcode scanner

- Automatictest recognitionvia2Dbarcodescanner
- Clear measurement process from cuvette insertion to storage
- Easy operation




The optics are insensitive to exter- nal light allowing measurements of MACHEREY-NAGEL cuvette tests without covering the cuvette slot. However, the cover is important for the performance of color measure- ments and measurements in the UV range.

Photometry becomes a pleasure

The all-rounders for all requirements

As complete spectrophotometers, the NANOCOLOR® VIS II and NANOCOLOR® UVIVIS II meet all requirements of your daily laboratory work. They come with the well-known barcode technology for NANOCOLOR® tube tests for fast measurements. Be- sides, they offer the possibility of extensive CIE-conform color measurements and scan-performance in real time. In addition to turbidity measurements, the devices offer comprehensive and free of charge brewery analysis according to MEBAK.

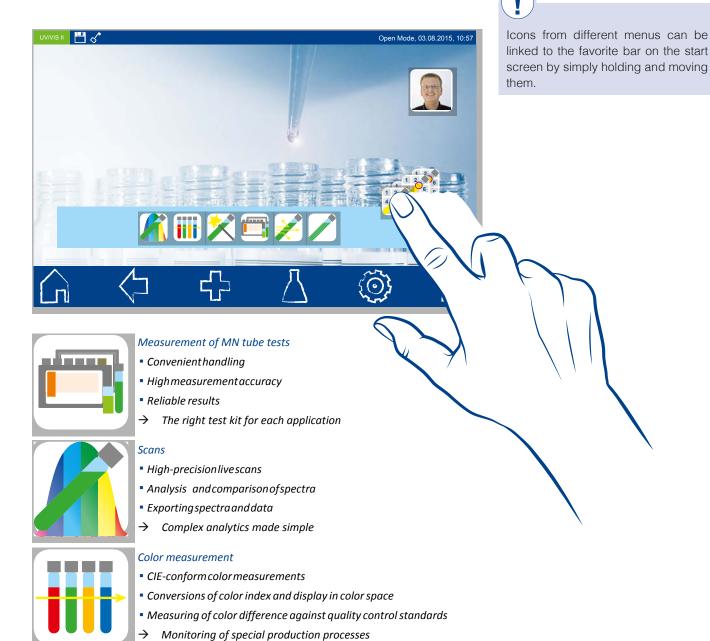




Turbidity measurement

- Turbidity measurement according to ENISO 7027
- Absorbance measurement (transmitted light measurement) at a 180 ° angle from 2–400 NTU
- Nephelometric stray light measurement from 0.1–1000 NTU
- → Fulfill requirements and assure measurement results

Photometry becomes a pleasure

The all-rounders for all requirements

NANOCOLOR® VIS II and NANOCOLOR® UV/vis II reduce the complexity of your daily laboratory work with easy menu guidance and icon-based procedures of the cuvette tests stored in the instrument. The clear result window enables an easy assignment of sample information to the measurement result. The step-by-step menu guidance for the calibration of special methods creates a simple entry level into user specific analytics even for inexperienced users.

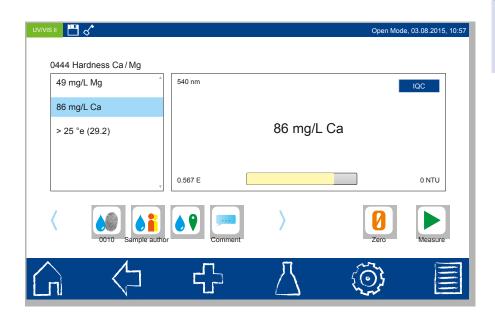



For each measurement up to eight sample information can be added and saved.



Result view

- Clearconnection of results and sample information
- Comprehensive presentation of all measurement results
- Direct linkage to IQC-options
- → All information at a glance



NTU-Check

- Safelyspot interfering turbidities
- Unique for maximum measurement result safety
- Runs automatically with each measurement of MN tube tests
- → Yourinsuranceforsecure measurement results




Pictograms in the instrument

- Easytounderstandpictograms
- Step-by-step procedures directly in the instrument
- Conduct measurements without extensive training
- → Reduce errors, enhance security



Method development

- Creation of users pecific special methods
- Comfortable calibration of customized methods
- Databaseforyourapplications
- → Versatilityforhighestdemands

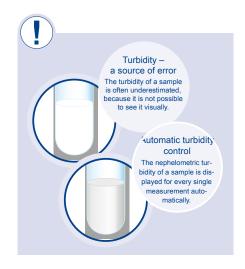

Be convincing

Your next audit will be a breeze

Quality is of high importance for MACHEREY-NAGEL. Thus, the brand-new spectro-photometers are equipped with extensive quality control features. Besides the integrated, f.o.c. inspection equipment monitoring, the devices offer a variety of quality control functions like standard measurements, multiple determinations and dilution series. In this way the NANOCOLOR® VIS II and NANOCOLOR® UV/vis II provide easy to use control options which enable an efficient and exact internal quality control.





Besides the monitoring of the complete analysis system the inspection equipment monitoring offers extensive possibilities for the inspection of instruments functionality. Perform the inspection equipment monitoring yourself and save money; additional external equipment checks are not necessary.



IQC

- IQC chart directly in the device
- Standard measurements, multiple determinations, dilution series and spike additions
- IQC memory for a clear documentation
- → Quality control at its best



User administration

- GLP-conform data storage
- Comprehensive user administration with password protection
- Individual adaptation to the respective user
- → Fulfill requirements



Inspection equipment monitoring

- Inspection equipment monitoring directly in the device
- Wavelengthcheckwith integratedholmiumoxidefilter
- Photometric accuracy check (NANOCONTROL NANOCHECK)
- → Ease controls and reduce costs

Stay connected

Striking interface options for smart connectivity

The connection of measuring devices to laboratory information systems (LMS) plays a more and more important role in many industries. Therefore, the NANOCOLOR® VIS II and NANOCOLOR® UV/vis II are equipped with all important interfaces (LAN, RS 232, USB) for the connection to modern laboratory information systems. Furthermore, the integrated LIMS configurator allows a customized adaptation for many kinds of data transfer.







LIMS configurator

- Easy implementation of measurement data to LIMS
- Customized composition of datasets
- → Tailor-madedataforindividualneeds



Interfaces (USB, LAN, RS 232)

- Future-proof interface options
- Straightforward data transfer to storage media
- Easy accessible USB-connections
- Future-proof connectivity options

Technical data

	NANOCOLOR ^{® UV} /vis II	NANOCOLOR® VIS II
Туре	Spectrophotometer with reference detector technology (RDT)	
Light sources	Halogen lamp (visible range) deuterium lamp (UV range)	Halogen lamp
Optical system	Monochromator	
Wavelength range	190–1100 nm	320–1100 nm
Wavelength accuracy	± 1 nm	
Wavelength resolution	0.1 nm	
Wavelength calibration	Automatic	
Wavelength selection	Automatic, barcode, manual	
Scan speed	1 complete scan in less than 1 min	
Spectral bandwidth	< 2 nm	< 4 nm
Photometric range	± 3.0 E in wavelength range 200–900 nm	± 3.0 E in wavelength range 340–900 nm
Photometric accuracy	0.005 E at 0.0-0.5 E; 1 % at 0.5-2.0 E	
Photometric linearity	< 0.5% at 2 E; ≤ 1% at > 2 E	
Stray light	< 0.05 %	< 0.1 %
Measuring modes	More than 200 preprogrammed tests, 100 optionally programmable methods, absorbance, transmission, factor, kinetics, 2-point calibration, scan, nephelometric turbidity measurement	
Turbidity measurement	Nephelometric turbidity measurement at 860 nm, 0.1–1000 NTU	
Cuvette slot	Test tubes 16 mm AD, rectangular cuvettes 2, 10, 20, 40, 50 mm	
Data memory	16 GB Micro SDHC card, 5000 measured data sets / spectra, GLP conform	
Display	10.1" HD display, anti-reflective cover glass with projective-capacitive touch screen (PCAP) and background LED lighting	
Operation	Barcode technology, icon-based menu guidance, touch screen	
Languages	DE/EN/FR/ES/PT/PL/HU	
External light	Insensitive, open cuvette slot	
Interfaces	LAN, 2 x USB (Host), 1 x USB (Function) and RS 232	
Update	Via Internet and USB-Stick	
Operating range	10–40 °C, max. 80 % relative humidity (without condensation)	
Power supply	110–240 V, ~50/60 Hz	12 V DC 3A
Dimensions L/W/H	400 / 440 / 170 mm	360 / 400 / 110 mm
Weight	6.5 kg	4.0 kg
Warranty	2 years	
Ordering information	Spectrophotometer NANOCOLOR® UV/vis II	Spectrophotometer NANOCOLOR® VIS II

REF 919 600

Spectrophotometer NANOCOLOR® V REF 919 650



incl. manual (Quickstart guide), protective covering, power cable, USB cable, USB stick, calibration cuvette, touchpen, cleaning cloth and certificate

incl. manual (Quickstart guide), protective covering, power supply, USB cable, USB stick, calibration cuvette, touchpen, cleaning cloth and certificate



Envitech Ltd

Unit S7 Capital Business Park Parkway Cardiff CF3 2PU Email: info@envitech.co.uk

Tel: 0044 2920 364252

Web: www.envitech.co.uk

